# Temperature dependent sensitivity of iodide chemical ionization mass spectrometers

Michael A. Robinson<sup>1,2,3</sup>, J. Andrew Neuman<sup>1,2</sup>, L. Gregory Huey<sup>4</sup>, James M. Roberts<sup>1</sup>, Steven S. Brown<sup>1,3</sup>, and Patrick R. Veres<sup>1</sup>

- <sup>5</sup> NOAA Chemical Sciences Laboratory, Boulder, Colorado, USA
  - <sup>2</sup>Cooperative Institute for Research in Environmental Sciences (CIRES), University of Colorado Boulder, Boulder, Colorado, USA
  - <sup>3</sup>Department of Chemistry, University of Colorado Boulder, Boulder, Colorado, USA
  - <sup>4</sup>School of Earth and Atmospheric Science, Georgia Institute of Technology, Atlanta, GA, USA
- 10 Correspondence to: Michael A. Robinson (michael.a.robinson@noaa.gov) and Patrick R. Veres (patrick.veres@noaa.gov)

#### **Abstract**

Iodide chemical ionization mass spectrometry (CIMS) is a common analytical tool used in both laboratory and field experiments to measure a large suite of atmospherically relevant compounds. Here, we describe a systematic ion molecule reactor (IMR) temperature dependence of iodide CIMS analyte sensitivity for a wide range of analytes in laboratory experiments. Weakly bound iodide clusters such as HCl, HONO, HCOOH, HCN, phenol, 2-nitrophenol and acyl peroxy nitrate (PAN) detected via the peroxy radical cluster, all exhibit strong IMR temperature dependence of sensitivity ranging from -3.4 to 5.9 % °C<sup>-1</sup> (from 37 to 47 °C). Strongly bound iodide clusters such as Br<sub>2</sub>, N<sub>2</sub>O<sub>5</sub>, ClNO<sub>2</sub>, and PAN detected via the carboxylate anion all exhibit little to no IMR temperature dependence ranging from 0.2 to -0.9 % °C<sup>-1</sup> (from 37 to 47 °C). The IMR temperature relationships of weakly bound clusters provide an estimate of net reaction enthalpy and comparison with database values, indicating these clusters are in thermal equilibrium. Ground site HCOOH data collected in the summer of 2021 in Pasadena, CA are corrected showing a reversal in the diel cycle, emphasizing the importance of this correction (35  $\pm$  6 % during the day, -26  $\pm$  2 % at night). Finally, we recommend two approaches to minimizing this effect in the field, namely heating or cooling the IMR; the latter has the added benefit of improving absolute sensitivity.

#### 1 Introduction

Atmospheric trace gases participate in chemical reactions that impact air quality and climate. Quantification of these trace gases, often at parts per trillion levels, demands highly sensitive and selective measurements. For decades, chemical ionization mass spectrometers (CIMS) have been used to provide rapid and sensitive in situ measurements of a wide variety of trace gases and aerosol composition. More recently, time of flight mass spectrometers (ToF-MS) have been developed to measure tens or even hundreds of compounds simultaneously, making their field deployment increasingly common in the last decade (Bertram et al., 2011; Lee et al., 2014; Lopez-Hilfiker et al., 2014; Krechmer et al., 2018; Koss et al., 2018; Lopez-Hilfiker et al., 2019; Veres et al., 2021). CIMS instruments employ reagent ions to react selectively with analytes of interest, and a number of

negative (I<sup>-</sup>, Br<sup>-</sup>, CF<sub>3</sub>O<sup>-</sup>, SF<sub>6</sub><sup>-</sup>, CH<sub>3</sub>C(O)O<sup>-</sup>) and positive (H<sub>3</sub>O<sup>+</sup>, NO<sup>+</sup>, NH<sub>4</sub><sup>+</sup>) reagent ions have been utilized in field and lab experiments to quantify a wide range of atmospheric trace gases (Slusher et al., 2004; Crounse et al., 2006; Huey et al., 1995; Veres et al., 2008; de Gouw and Warneke, 2007; Koss et al., 2016; Zaytsev et al., 2019). Negative ion-adduct chemistry has been used with quadrupole mass spectrometers that sample the atmosphere since 1975 (Dougherty et al., 1975). The use of iodide-adduct ion chemistry has become commonplace in atmospheric chemistry field campaigns as this chemistry has provided sensitive measurements of molecular halogens, halogen oxides, carboxylic acids, acyl peroxynitrates (PANs), oxygenated organics and nitrogen oxides (Slusher et al., 2004; Osthoff et al., 2008; Kercher et al., 2009; Thornton et al., 2010; Mielke et al., 2011; Riedel et al., 2012; Neuman et al., 2016).

Often iodide ions (and/or its water clusters) react with a neutral analyte molecule to form an adduct that can be identified and detected in a ToF mass spectrometer. Iodide's large mass defect, low likelihood of electron transfer or proton abstraction often allow selective ionization chemistry that can achieve low detection limits and avoid interferences. Additionally, iodide's single natural isotope greatly simplifies complex mass spectral interpretation. Iodide adduct formation is typically only slightly exergonic (i.e. ΔG < 0), meaning adduct formation is rapid (i.e. every collision forms a cluster ion), and relatively fragment free (i.e. soft ionization) (de Hoffmann and Stroobant, 2007). A complication of iodide CIMS, in contrast to proton transfer reactions, is that sensitivities to different analytes range several orders of magnitude (Iyer et al., 2016) which implies that product ion abundance is not solely determined by forward reaction rates. It has been proposed that this range of sensitivities is due to the variation in binding enthalpy between the iodide anion and neutrals (Iyer et al., 2016).

The iodide anion-molecule chemistry in the presence of water is a system of reactions:

$$I^- + A \rightleftharpoons [I \cdot A]^- \tag{R1}$$

$$[I \cdot H_2 O]^- + A \rightleftharpoons [I \cdot A]^- + H_2 O$$
 (R2)

$$I^- + RC(0)O_2 \to RC(0)O^- + IO$$
 (R3)

Here, I<sup>-</sup> and I<sup>-</sup> •H<sub>2</sub>O are reagent ions, A is the analyte of interest, and I•A<sup>-</sup> is a product ion. Reaction R3 pertains to the special analyte case of peroxy acetyl radicals, which differs from most analyte cluster reaction pathways (i.e. R1 or R2). The portion of each reaction pathway depends on the analyte and reagent ([I•H<sub>2</sub>O<sup>-</sup>] or [I<sup>-</sup>]). Ambient mixing ratios of a given analyte are obtained using Eq. 1:

$$[A] = \frac{[I \cdot A^{-}]}{[I^{-}]kt} \tag{1}$$

Where k is the effective rate constant of the reaction and t is the reaction time, and this equation also holds true for other reagent ions (i.e.  $I^- \bullet H_2O$ ). Typically k and t are not determined separately, but the product of the two parameters, or the effective sensitivity, is determined for a given analyte (Huey, 2007).

Much has been done to characterize this ion chemistry for application to atmospheric chemistry field and laboratory studies (Huey et al., 1995; Slusher et al., 2004; Kercher et al., 2009; Lee et al., 2014; Lopez-Hilfiker et al., 2016). Analyte ionization occurs subsequent to mixing with reagent ions in a volume known as the ion molecule reactor (IMR), see Figure 1. Instrument analyte detection sensitivity is impacted by the effective reaction rate constants and pathways of the analyte and reagent ion(s)

in the IMR (R1-R3). Water vapor plays an important role in the iodide-adduct ionization efficiency, which in turn impacts instrument sensitivity (Lee et al., 2014; Slusher et al., 2004; Huey et al., 1995; Kercher et al., 2009). Instrument sensitivity for a given analyte molecule can exhibit a strong to no relationship with water vapor, and this water dependence requires careful calibration of these instruments as a function of IMR water vapor pressure (Lee et al., 2014). This instrument sensitivity dependence on water has motivated efforts to actively control water or carefully calibrate as a function of water in order to stabilize or determine sensitivity during field deployment (Lee et al., 2014; Slusher et al., 2004; Veres et al., 2021). IMR temperature determines the importance of kinetics versus thermodynamics for determining the outcome of ion molecule reactions, yet there is limited discussion of the importance of IMR temperature in CIMS systems.

Reactions in the IMR can be in a kinetic or thermodynamic control regime, where the kinetic regime is described as a sufficiently exergonic reaction (i.e.  $\Delta G < 0$ ) with a slow reverse reaction rate, such that equilibrium is not established in the IMR. The thermodynamic control regime is described as moderately exergonic reaction (i.e.  $\Delta G < 0$ ) with reverse reaction rates ~5 times smaller than the residence time in the IMR allowing the ion-analyte system to reach equilibrium on this timescale. The SiF<sub>5</sub>-(HNO<sub>3</sub>) system achieves equilibrium on timescales equivalent to the residence time in most IMRs used in field deployable CIMS (~100 ms) (Huey, 2007). A CIMS using Br reagent ions to measure the Br (HO<sub>2</sub>) cluster exhibited thermal equilibrium behavior, with a decrease in instrument sensitivity of 20% with an increase of IMR temperature from 20 to 40 °C (Sanchez et al., 2016). Similar thermochemistry systems (i.e. similar  $\Delta G$ ) are likely to equilibrate on similar timescales.

Limited kinetic studies have been conducted for iodide anion – molecule reactions. Examples include the  $N_2O_5$  clustering reaction (R1) (k = 1.3 x  $10^{-9}$  cm<sup>3</sup> s<sup>-1</sup> at 293 K) and the carboxylate anion reaction (R3) (k = 4 x  $10^{-10}$  to 16 x  $10^{-10}$  cm<sup>3</sup> s<sup>-1</sup> at 298 K) (Huey et al., 1995, p.6; Villalta and Howard, 1996). However, quantum chemical calculations predict collision rates of neutral molecules with I<sup>-</sup> for a variety of atmospherically relevant molecules to range from 2 – 4 x  $10^{-9}$  cm<sup>3</sup> s<sup>-1</sup> (Iyer et al., 2016). If an analyte is in the thermal controlled regime, the temperature dependent equilibrium constant of this system of reactions will affect the instrument sensitivity. A system in equilibrium is described by the linear form of the Van't Hoff equation which relates the change in equilibrium constant ( $K_{eq}$ ) of the ion-molecule reaction to the absolute temperature (T):

$$\ln K_{eq} = -\frac{\Delta_r H^0}{RT} + \frac{\Delta_r S^0}{R} \tag{2}$$

where  $\Delta_r H^0$  is the standard reaction enthalpy,  $\Delta_r S^0$  is the standard reaction entropy, and R is the gas constant. The slope and intercept of the Van't Hoff plot (ln K<sub>eq</sub> vs 1/T) gives the enthalpy and entropy change, respectively (Keesee and Castleman, Jr., 1986).

Analyte sensitivity in CIMS is a product of the net rate of product ion formation, a function of adduct binding energy, and the instrument transmission of those adduct ions. A recent framework for understanding iodide adduct binding energies uses a mass spectrometer ion lens voltage scanning procedure to estimate the relative binding energies (Lopez-Hilfiker et al., 2016). This approach gives important insight to the stability of ion adducts and the extent to which they may decluster during transmission, and can be used to provide inferred or relative sensitivities (Lopez-Hilfiker et al., 2016). A drawback of the

95

voltage scanning method for determination of relative adduct binding energies is the potential for unintended changes in instrument ion transmission, which may result from modulation of instrument tuning to induce the needed change in field strength. In this work we leverage kinetic energy in the form of a temperature change as a substitute for modulation of the field strength, thereby eliminating the impact of ion transmission. The effect of temperature on instrument sensitivity has not been widely explored for the iodide adduct CIMS, where the impact has been documented for limited analytes such as the iodide carboxylate anion reaction (R3) and the  $Br^-(HO_2)$  cluster (Villalta and Howard, 1996; Sanchez et al., 2016).

Previously, several CIMS IMR temperature control methods have been implemented. Documented instrument deployments range from actively heating the IMR region or hardware close to the IMR (50 to 60 °C), such as for the Filter Inlet for Gases and Aerosols instrument (Lopez-Hilfiker et al., 2014), to insulation of the IMR region (Lee et al., 2018), to active cooling of the IMR region (15 °C) (Neuman et al., 2002). However, thermal coupling between the IMR and ToF body, which may be difficult to temperature control in field deployments of CIMS instruments (discussed below), may introduce temperature variations in regions where ion molecule reactions occur.

110 We extend the impact of temperature on the abundance of product ions to a range of atmospheric trace gases and provide a framework for understanding variations of detection sensitivity with temperature for I<sup>-</sup> ion molecule reactions. Additionally we estimate thermally controlled reaction pathways in our IMR system by comparing observed reaction enthalpies to literature values. Finally, we compare our experiments to field observations and recommend temperature control strategies to improve sensitivity stability.

#### 115 **2** Experimental Approach

100

105

120

The experiments here are conducted on a modified commercially available CIMS (https://www.tofwerk.com/products/api-tof/) that has been used extensively for sampling from ground and airborne platforms (Veres et al., 2021). Figure 1 depicts the front end of the instrument, which differs from the commercially available instrument. The NOAA I<sup>-</sup> ToF-CIMS (hereafter NOAA CIMS) utilizes a custom pressure controlled inlet that maintains constant instrument pressures independent of ambient pressure. The NOAA CIMS notably operates at lower IMR pressure (40 mbar, rather than 100 mbar used in many instruments) and small segmented quadrupole (SSQ) pressure (1.64 mbar, rather than 2 mbar) in order to reduce the impact of secondary chemistry in the IMR and SSQ (Bertram et al., 2011; Lee et al., 2014). The remainder of the ToF has not been modified.

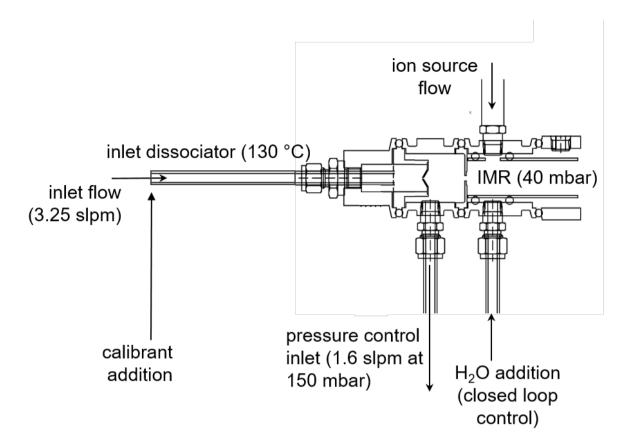


Figure 1: Instrument inlet schematic.

125 The NOAA CIMS detects analyte molecules with I<sup>-</sup> and I<sup>-</sup>•H<sub>2</sub>O ion clusters via ToF-MS with a mass resolving power of approximately 5000 (Δm/m). Reagent ions are formed by flowing 1 standard liter per minute (slpm) of N<sub>2</sub> and 9 standard cubic centimeter per minute (sccm) of 1025 ppm methyl iodide (CH<sub>3</sub>I in N<sub>2</sub>) in front of a vacuum ultraviolet Krypton lamp (Restek photoionization lamp, model 108-BTEX, Restek Corporation, PA, USA). This ionization technique produces an intense and relatively interference-free source of iodide ions (Ji et al., 2020; Breitenlechner et al., 2022). These experiments (unless stated otherwise) were done with water dynamically added to the IMR (via a saturated small N<sub>2</sub> flow (0 – 100 sccm)) to maintain a cluster ratio (I<sup>-</sup>•H<sub>2</sub>O:I<sup>-</sup>) (CR) of 0.50 ± 0.002, as is typical in field operation of the instrument. The added saturated N<sub>2</sub> flow can dilute the sample in the IMR and is typically a negligible effect, but due to the nature of these experiments (large N<sub>2</sub> flow at high IMR temperatures) we have corrected for dilution. A detailed description of the CR closed loop control system can be found in the Supplement Section 1.

135 Ambient air was sampled through a mass flow controlled (3.25 slpm) PFA inlet (30 cm length, 0.64 cm ID), which was heated to 130 °C to thermally dissociate PAN for detection via the carboxylate anion (R3) or peroxy radical cluster (Slusher et al., 2004). Ambient air is sampled through a 1000 μm diameter orifice into a pressure controlled region maintained at 150 mbar. This pressure control guarantees that 1.2 slpm enters the IMR, independent of ambient pressure, and mixes with the ~1 slpm

flow from the ion source. The reagent ion signals during this work were typically 1×10<sup>6</sup> ion counts per second (MHz) for I•H<sub>2</sub>O and 2 MHz for I<sup>-</sup>. Ion packet extraction frequency in the time of flight region of the instrument during this work was
25,000 Hz. Product cluster ions are normalized by the I<sup>-</sup>•H<sub>2</sub>O signal to account for changes in ion source output. Every two
minutes, the calibrants are removed and lab air sampled, and the calibration signals are determined as the difference between
lab air alone and lab air with calibrant. Two IMR designs are investigated in this work, the standard Aerodyne Inc, IMR (ARI)
and the NOAA IMR, which has been deployed on the instrument for the Atmospheric Tomography – 3 aircraft campaign
(ATOM – 3, 2017) (https://espo.nasa.gov/atom/archive/browse/atom/id10), ATOM – 4 (aircraft campaign, 2018)
(https://espo.nasa.gov/atom/archive/browse/atom/id14), Fire Influence on Regional to Global Environments and Air Quality
(FIREX-AQ) (aircraft campaign, 2019) (https://www-air.larc.nasa.gov/missions/firex-aq/) and Southwest Urban NOx and
VOC experiment (SUNVEx) (ground site, 2021) (https://csl.noaa.gov/projects/sunvex/). The NOAA IMR design was
deployed during SUNVEx and the ARI IMR design was used to investigate the influence of IMR temperature in the laboratory.
Both IMR designs are summarized in Supplement Table S1.

#### 2.1 IMR Temperature experiments

140

145

150

155

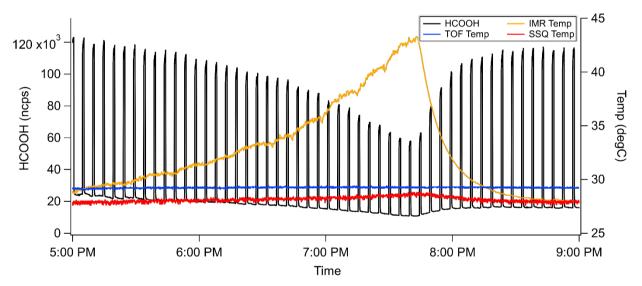


Figure 2: Typical time series for an IMR temperature experiment. Product ion signal is normalized to I-•H<sub>2</sub>O signal.

The ion chemistry that forms unique product ions occurs primarily in the IMR. Consequently, the temperature dependence of the sensitivity was predominately examined as a function of IMR temperature. Experiments were conducted using the widely used temperature controllable ARI IMR, which includes a cartridge heater in a stainless steel block mounted to the IMR and controlled via the instrument computer system (Bertram et al., 2011; Krechmer et al., 2016). All temperatures are obtained from manufacturer supplied thermistors placed on the IMR and ToF. Even without independently heating the IMR, the heat generated by the instrument maintains the IMR temperature at approximately 30 °C. IMR temperature was increased from 30

°C to 45°C over a period of three hours as shown in Figure 2. After reaching the maximum IMR temperature, temperature control was turned off to allow the IMR to return to ambient temperature. During the entire experiment, standard additions of analyte molecules occurred every 2 minutes for 45 seconds. Normalized product ion signal and instrument temperatures during a typical IMR temperature experiment with formic acid (HCOOH) additions is shown in Figure 2.

#### 2.2 ToF Body temperature experiments

165 As discussed in section 3.3 below, temperature dependent ion chemistry can occur downstream from the IMR, in which case IMR temperature alone does not determine sensitivity. These regions contain segmented quadrupoles that guide ions, and their temperature is determined largely by the temperature of the entire ToF body, which is closely coupled with ambient temperature. Collisional dissociation occurs during voltage scanning experiments which vary voltages between the first (SSQ) and second (big segmented quadrupole (BSQ)) quadrupole ion guides, indicating SSQ temperatures may affect sensitivity 170 (Lopez-Hilfiker et al., 2016). While newer instruments have improved temperature control of low pressure regions of the instrument (IMR, SSQ, BSQ and ToF), for many applications it may not be feasible to control the temperature of the entire ToF, since it is a very large thermal mass that is often securely mounted when installed on a mobile platform. Consequently, we have investigated the impact of a variable ToF temperature, while controlling IMR temperature, to represent typical field conditions. These additional experiments were conducted using a radiative heater placed next to the ToF body and insulating 175 blanket. Three types of experiments were performed using the ARI IMR design to examine temperature control strategies that could be implemented to stabilize sensitivity if ambient temperatures change substantially. The first is to actively cool the IMR and determine sensitivity while the ambient temperature was varied. This control strategy may afford improved sensitivity for many molecules and reduces sensitivity changes caused by large changes in ambient temperature. The IMR was cooled using a thermoelectric chiller that pumps cold water through copper tubing wrapped around the IMR. This control strategy achieved 180 an IMR temperature of ~18 °C during these experiments. The second control strategy is to heat the IMR to a temperature higher than the maximum expected ToF body temperature (in this case ~45 °C). This control strategy suffers from a loss in sensitivity for some analytes, but is more easily implemented using the commercially provided hardware compared to cooling the IMR. For both of these active control strategies, the IMR temperature was allowed to stabilize, then the radiative heater was turned on to a set point of ~40 °C, which replicates many harsh field deployment temperatures in aircraft and at ground 185 sites. The third strategy is the absence of active temperature control. This experiment replicates conditions during prior field deployments and reveals the sensitivity changes that can occur in the absence of IMR temperature control and was conducted on both IMR designs. Product ion signals, reagent ion signals and temperatures during a typical experiment without IMR temperature control are shown in Supplement Figure S2. A comparison of IMR saturated N<sub>2</sub> flow rates during the ToF body temperature experiment between the two investigated IMR designs is provided in Supplement Figure S3. The changes observed 190 in added IMR saturated N<sub>2</sub> flow is likely due to geometry differences between the two designs and the impact on mixing in the IMR.

## 2.3 Calibration sources

195

200

In order to understand the broader impacts of IMR temperature on detection sensitivity, a wide range of analytes (see Table 1) were added to the instrument while varying the IMR temperature. Iodide ion chemistry detection sensitivity and water dependence varies considerably among analytes (for example Lee et al, 2014, SI). We have tested the IMR temperature dependence for compounds that span a wide range of sensitivities. Table 2 summarizes the calibration method, inlet mixing ratio, reagent-normalized instrument sensitivity (CR = 0.5, IMR temperature = 29 °C) and sensitivity dependence on IMR temperature. Permeation tubes (HCOOH and Br<sub>2</sub>: Kin-Tek; HCl: VICI Metronics) are held in PFA sleeves and temperature controlled (HCl and Br<sub>2</sub> at 30 °C, HCOOH at 40 °C). A continuous 50 sccm flow of N<sub>2</sub> is directed around the permeation tubes and through the PFA sleeves. Phenol (99% MilliporeSigma) and 2-Nitrophenol (98% Sigma Aldrich) were added to a glass diffusion cell (Allen Scientific Glass) with a headspace flow of approximately 30 sccm zero air.

Analyte	Calibration Source	Inlet	Instrument Sensitivity (Hz pptv <sup>-1</sup>	Sensitivity
		Mixing	1 MHz <sup>-1</sup> I <sup>-</sup> •H <sub>2</sub> O) at CR = 0.50 and	dependence
		ratio	IMR temperature = 29 °C	on IMR
				Temperature
				(% °C <sup>-1</sup> ) <sup>a</sup>
HCN	10.4 ppmv in N <sub>2</sub> gas	20 ppbv	0.05	$-5.7 \pm 0.3$
	standard			
HCl	Permeation tube	27 ppbv	0.0013	$-5.8 \pm 0.4$ 205
НСООН	Permeation tube	7 ppbv	5.3	$-5.0 \pm 0.2$
$Br_2$	Permeation tube	1.3 ppbv	10.8	$-0.79 \pm 0.05$
ClNO <sub>2</sub>	Online synthesis	30 pptv	11	$-0.89 \pm 0.4$
N <sub>2</sub> O <sub>5</sub>	Online synthesis	3.8 ppbv	36	$-0.90 \pm 0.1$
HONO	Online synthesis	2.5 ppbv	2.3	$-4.47 \pm 0.6$
I*PAN	Online synthesis	800 pptv	2.4	$-3.41 \pm 0.1$
PAN Anion	Online synthesis	800 pptv	2.1	$0.2 \pm 0.09$
2-Nitrophenol	Diffusion cell	350 pptv	0.01 (Palm et al., 2020)	$-5.9 \pm 0.3$
Phenol	Diffusion cell	1055 pptv	0.58 (Palm et al., 2020)	$-5.4 \pm 0.4$

210

Table 1: Calibration sources, instrument sensitivity and IMR temperature dependence of sensitivity for analytes studied.

<sup>&</sup>lt;sup>a</sup>Change in sensitivity from 310 to 320 K

## 3 Results and Discussion

215

220

225

230

Most analytes detected in the NOAA CIMS lose sensitivity with increasing temperature, shown in Figure 2 and Table 2. Notably the sensitivity of several halogens remain nearly constant with temperature and the PAN carboxylate anion (R3) sensitivity does not decrease with temperature. Instrument sensitivities are normalized to the sensitivity measured at 310 K, as this allows for a comparison across all field and lab operating conditions examined here. Additionally, normalization permits a comparison of the temperature dependence of analytes with absolute sensitivities (Hz pptv<sup>-1</sup>) that range several orders of magnitude between analytes (Lee et al., 2014; Iyer et al., 2016).

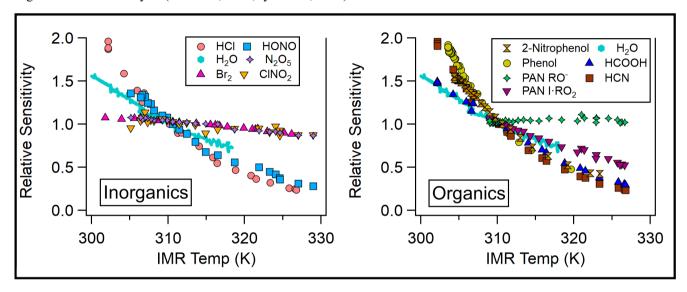


Figure 2: Relative sensitivity (normalized to sensitivity at 310 K) as a function of IMR temperature for organic and inorganic atmospheric trace gases detected via Iodide adduct chemical ionization.

The analytes studied here can be grouped into two classes, strongly bound clusters with weak IMR temperature dependence (such as Br<sub>2</sub>, N<sub>2</sub>O<sub>5</sub>, and ClNO<sub>2</sub>) and weakly bound clusters with strong IMR temperature dependence (HCl, HONO, HCOOH, HCN and PAN via peroxy radical cluster). The carboxylate anion detection of PAN behaves similarly to strongly bound clusters, but is not clustering ion-molecule chemistry and cannot be classified. The absolute sensitivities of Phenol and 2-Nitrophenol are currently unknown for our instrument, so we have not classified them. However we expect them to have low sensitivity, since these molecules are detected with relatively low sensitivity in other iodide CIMS (0.58 and 0.01 Hz pptv<sup>-1</sup> respectively) (Palm et al., 2020). Weakly bound inorganic analytes, such as HCl and HONO exhibit reductions of sensitivity of 80% and 67% respectively for an IMR temperature increase from 31 to 50 °C. Weakly bound organic analytes, such as Phenol and HCOOH, exhibit decreases in sensitivity of 75% and 60% respectively when the IMR temperature increases from 30 to 46 °C. Strongly bound inorganic analytes, such as Br<sub>2</sub> and N<sub>2</sub>O<sub>5</sub>, exhibit much smaller changes in sensitivity (15%) over these IMR temperature ranges.

#### 3.1 Impact of IMR heating on water cluster distribution

IMR heating impacts both product ion formation and the iodide water cluster distribution. We believe that iodide and its first water cluster (I<sup>-</sup>•H<sub>2</sub>O) dominate the reagent ion distribution in our IMR. The detector signal for the most abundant higher order cluster is the second water cluster (I<sup>-</sup>•(H<sub>2</sub>O)<sub>2</sub>), which is typically 0.05% of I<sup>-</sup>•H<sub>2</sub>O, and likely reflects the low abundance of higher order clusters in the IMR. IMR temperature sweep experiments were conducted to an IMR temperature of ~55 °C. This is the maximum temperature in which instrument CR could be controlled to 50 % which is the nominal operating CR during field deployment. The reduction in the iodide water cluster product ion as a function of IMR temperature is shown in Figure 2. This experiment was conducted without IMR temperature or water control, in order to understand the influence of IMR temperature on the I<sup>-</sup>•(H<sub>2</sub>O) product ion. This thermodynamically controlled behavior of the I<sup>-</sup>•(H<sub>2</sub>O) product ion, which is also a critical reagent ion in our IMR, emphasizes the importance of temperature and CR control to maintain a constant reagent ion distribution and stable sensitivity.

#### 3.2 IMR temperature dependence gives insight into iodide ion chemistry

245

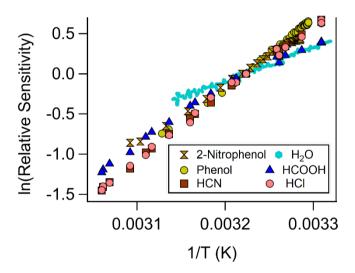


Figure 3: Van't Hoff relationship for a selection of atmospheric trace gases using Iodide adduct ion chemistry in the NOAA CIMS (normalized at 310 K).

Analyte systems that equilibrate in the IMR can be investigated further with the Van't Hoff relationship (Keesee and Castleman, Jr., 1986), shown in Figure 3. The positive slopes indicate that the net reactions of these analytes are exothermic (i.e. Δ<sub>r</sub>H < 0). Analyte systems in equilibrium give a measurement of the net reaction enthalpy in the IMR, according to the slope of the Van't Hoff relationship (Eq. 2) in Figure 3 (see Supplement Figure S4 for HCN example). These measured reaction enthalpies can be represented by a combination of the two predominant reaction pathways (R1 & R2) reaction enthalpies. These condition-specific reaction enthalpies are compared to reported values (NIST Chemistry WebBook, SRD 69, 2021),

shown in Figure 4. In order to compare measured reaction enthalpies to literature values, the ligand switching (R2) reaction enthalpies must be calculated, by subtracting the single water cluster reaction enthalpy from the analyte direct cluster reaction enthalpy (blue in Figure 4), an example calculation is shown below for the HCOOH ligand switching reaction:

260

$$[H_2O \cdot I^-] \rightleftharpoons H_2O + I^-$$

$$HCOOH + I^- \rightleftharpoons [HCOOH \cdot I^-]$$

$$LCOOH + [H_2O \cdot I^-] \rightleftharpoons [HCOOH \cdot I^-] + H_2O$$

$$\Delta_r H^0 = -79.1 \pm 4.2 (kJ/mol)$$

$$\Delta_r H^0 = -79.1 \pm 4.2 + 43 \pm 3 = -36.1 \pm 7.2 (kJ/mol)$$

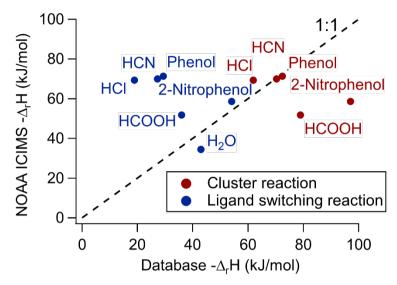


Figure 4: NOAA CIMS measured reaction enthalpies vs NIST webbook database reaction enthalpies for either the direct clustering pathway or ligand switching pathway. A 1:1 line is included to indicate where the measured reaction enthalpy agrees with the database.

265 The direct cluster and ligand switching reaction enthalpies are shown in red and blue, respectively in Figure 4. Where the 1:1 line falls between these values is an indication of the apparent branching ratio between each reaction pathway in the IMR under these conditions. For the case of HCOOH, the net reaction enthalpy is  $-51.8 \pm 1.4$  kJ/mol. Consequently, this thermochemical system lies between the ligand switching reaction enthalpy of -36.1 ± 5.1 kJ/mol and the direct cluster reaction enthalpy of - $79.1 \pm 4.2$  kJ/mol. This result, which is only valid for conditions in the IMR during this experiment, implies that HCOOH 270 reacts both through the ligand switching pathway and the direct cluster pathway. A similar conclusion can be drawn for 2nitrophenol, as the measured net reaction enthalpy of  $-58.6 \pm 0.6$  kJ/mol falls between the reported ligand switching reaction enthalpy of -54.1  $\pm$  8.1 kJ/mol and direct clustering reaction enthalpy of -97.1  $\pm$  7.5 kJ/mol. These analytes differ from HCN, HCl and phenol which, within experimental uncertainty, all fall close to the 1:1 line for the direct clustering pathway (see Figure 4). For example, the literature reported value for the HCl clustering reaction is  $-61.9 \pm 8.4$  kJ/mol which compares well 275 within the experimentally determined value of  $-69.4 \pm 0.9$  kJ/mol. The involvement of the ligand switching chemistry may explain the relatively high sensitivity for HCOOH as compared to extremely low sensitivity direct clustering analytes, such as HCN, HCl and phenol.

## 3.3 Field deployment and IMR temperature control

280

285

290

295

The impact on HCOOH sensitivity due to IMR temperature changes during the SUNVEx field campaign (NOAA IMR) are shown in Figure 5. During the day (09:00 to 18:00 Pacific Daylight Time (PDT)), elevated instrument temperatures caused a decrease in HCOOH sensitivity, requiring campaign average correction of  $35 \pm 6$  % (average  $\pm 1$  standard deviation). At night (20:00 to 06:00 PDT), instrument temperatures were lower than lab determined sensitivities increasing HCOOH sensitivity, requiring campaign average correction of  $-26 \pm 2$  %. The diel dependence of HCOOH is reversed if this temperature correction is not applied making it critical to accurate interpretation of these data.

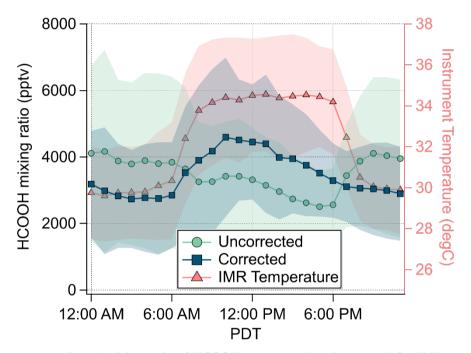


Figure 5: Campaign average diurnal mixing ratio of HCOOH (uncorrected and corrected for IMR temperature) and IMR temperature. Shading represents 1 standard deviation.

Sufficient energy can de-cluster even the most strongly bound clusters, as shown by this work and SSQ voltage scanning experiments (Lopez-Hilfiker et al., 2016). Therefore, in order to maintain stable sensitivities in the instrument, understanding the locations where declustering occurs is important to implementing efficient IMR temperature control strategies in the field. During the SUNVEx field deployment of the NOAA CIMS in 2021, the instrument was installed at a ground site, where the IMR temperatures ranged from 25 to 40 °C (298 to 313 K) (see Supplement Figure S5). The temperature changes were caused both by the PAN inlet dissociator that was on only during the day and ambient temperature fluctuations. This temperature range is not atypical for field locations (see Supplement Figure S6 for research aircraft example). Hourly HCOOH standard additions indicated a significant change in instrument sensitivity as a function of IMR temperature, shown in Figure 6b. The variations in sensitivity with temperature in lab experiments and field observations are similar, however absolute sensitivities differ between the two IMR designs (SUNVEx, NOAA IMR, shown in gray in Fig 6b; Lab experiments, ARI IMR, red, salmon

and blue points, Fig 6b). Small differences in these slopes indicate IMR temperature laboratory experiments do not fully mimic the behavior of the instrument when deployed in a harsh environment, suggesting that ion-analyte de-clustering occurs not only in the IMR, but further into the instrument, likely in the SSO or BSO.

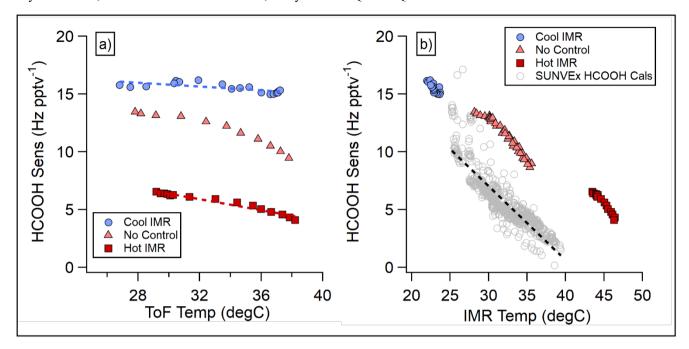


Figure 6: Comparison of absolute HCOOH sensitivity between laboratory IMR temperature experiments (ARI IMR design) with SUNVEx field standard additions (NOAA IMR design).

Several ToF body experiments were conducted on the ARI IMR to determine effective methods for stabilizing sensitivity when environmental temperatures change (Figure 6). Instrument temperatures (described as ToF temperature) varied from 28 to 38 °C over the course of these experiments, as shown in panel 6a. The heated IMR control strategy attempted to maintain the IMR at 45 °C while the ToF temperature varied. The cooled IMR control strategy attempted to maintain the IMR temperature at 18 °C while the ToF temperature varied. Both IMR temperature control strategies improve HCOOH sensitivity stability compared to an uncontrolled instrument. Although, HCOOH sensitivity is reduced with the heated IMR control (red squares), this control strategy still reduces sensitivity drift (~50% reduction) when compared to an uncontrolled IMR (salmon triangles). Cooling the IMR (blue circles) both improves absolute sensitivity and reduces sensitivity drift further (~75% reduction in sensitivity drift) when compared to uncontrolled IMR temperature operation.

## 4 Conclusions

300

305

310

315

Iodide ion chemistry, commonly used in ToF-CIMS that sample the atmosphere, imparts a temperature dependence to analyte sensitivity. The sensitivity variation with temperature for some analytes is large (5% °C<sup>-1</sup> from 36 to 40 °C). Importantly sensitivity can be largely controlled by temperature stabilizing the IMR. Further, the dependence of sensitivity on temperature

gives insights into the rates of ion molecule reactions that control the product ion abundance, and provides information on absolute sensitivity. Measuring the temperature dependence of the sensitivity allows for the determination of iodide analyte reaction enthalpies. These enthalpies reveal the reactions that control the abundance of product ions in the IMR. Weakly bound clusters are strongly temperature dependent and can be detected more sensitively at lower IMR temperatures. This work was prompted by the regular addition of a weakly bound analyte (HCOOH), which has a sensitivity temperature dependence similar to other weakly bound analytes. We recommend a similar regular calibration approach is taken for all CIMS field deployments. These considerations are likely not unique to the iodide ion adduct system but are more generally applicable to other CIMS ion chemistries. We recommend careful examination of IMR temperature dependencies of instrument sensitivities, especially for weakly bound clusters (many organics) where the effect may be large. Temperature control of the IMR region can help reduce the impact, but further de-clustering (either due to temperature or collisional fragmentation) of the product ion may occur in the instrument ion optics. It is worth noting that higher pressure IMR systems may be more susceptible to thermal effects due to the increased time to reach equilibrium. Future work on expanding the range of analyte classes (for example highly functionalized organics) and ion chemistries will be of interest to improve our skill at predicting CIMS absolute sensitivity.

## Data availability

320

325

330

340

Data from the SUNVEx campaign are available to the general public at <a href="https://csl.noaa.gov/groups/csl7/measurements/2021sunvex/GroundLA/DataDownload/">https://csl.noaa.gov/groups/csl7/measurements/2021sunvex/GroundLA/DataDownload/</a>

## **Author contribution**

MR, AN and PV designed the experiments, MR and AN carried them out. MR, AN and PV executed the field project. MR and AN prepared the manuscript with contributions from all co-authors.

### Acknowledgements

The NOAA Chemical Sciences Laboratory acknowledges support for this work from the California Air Resources Board under agreement number 20RD002. This work was supported in part by the NOAA Cooperative Agreement with CIRES, NA17OAR4320101.

## References

- Bertram, T. H., Kimmel, J. R., Crisp, T. A., Ryder, O. S., Yatavelli, R. L. N., Thornton, J. A., Cubison, M. J., Gonin, M., and Worsnop, D. R.: A field-deployable, chemical ionization time-of-flight mass spectrometer, Atmos. Meas. Tech., 4, 1471–1479, https://doi.org/10.5194/amt-4-1471-2011, 2011.
- 345 Breitenlechner, M., Novak, G. A., Neuman, J. A., Rollins, A. W., and Veres, P. R.: A versatile vacuum ultraviolet ion source for reduced pressure bipolar chemical ionization mass spectrometry, Atmos. Meas. Tech., 15, 1159–1169, https://doi.org/10.5194/amt-15-1159-2022, 2022.
  - Crounse, J. D., McKinney, K. A., Kwan, A. J., and Wennberg, P. O.: Measurement of Gas-Phase Hydroperoxides by Chemical Ionization Mass Spectrometry, Anal. Chem., 78, 6726–6732, https://doi.org/10.1021/ac0604235, 2006.
- Dougherty, R. C., Roberts, J. D., and Biros, F. J.: Positive and negative chemical ionization mass spectra of some aromatic chlorinated pesticides, Anal. Chem., 47, 54–59, https://doi.org/10.1021/ac60351a054, 1975.
  - de Gouw, J. and Warneke, C.: Measurements of volatile organic compounds in the earth's atmosphere using proton-transfer-reaction mass spectrometry, Mass Spectrom. Rev., 26, 223–257, https://doi.org/10.1002/mas.20119, 2007.
  - de Hoffmann, E. and Stroobant, V.: Ion Sources: Gas-Phase Ion-Molecule Reactions, in: Mass Spectrometry: Principles and
- 355 Applications, John Wiley & Sons, 72–76, 2007.
  - Huey, G. L., Hanson, D. R., and Howard, C. J.: Reactions of SF6- and I- with Atmospheric Trace Gases, J. Phys. Chem., 99, 5001–5008, https://doi.org/10.1021/j100014a021, 1995.
  - Huey, L. G.: Measurement of trace atmospheric species by chemical ionization mass spectrometry: Speciation of reactive nitrogen and future directions, Mass Spectrom. Rev., 26, 166–184, https://doi.org/10.1002/mas.20118, 2007.
- 360 Iyer, S., Lopez-Hilfiker, F., Lee, B. H., Thornton, J. A., and Kurtén, T.: Modeling the Detection of Organic and Inorganic Compounds Using Iodide-Based Chemical Ionization, J. Phys. Chem. A, 120, 576–587, https://doi.org/10.1021/acs.jpca.5b09837, 2016.
  - Ji, Y., Huey, L. G., Tanner, D. J., Lee, Y. R., Veres, P. R., Neuman, J. A., Wang, Y., and Wang, X.: A vacuum ultraviolet ion source (VUV-IS) for iodide–chemical ionization mass spectrometry: a substitute for radioactive ion sources, Atmos. Meas.
- 365 Tech., 13, 3683–3696, https://doi.org/10.5194/amt-13-3683-2020, 2020.
  - Keesee, R. G. and Castleman, Jr., A. W.: Thermochemical Data on Gas-Phase Ion-Molecule Association and Clustering Reactions, J. Phys. Chem. Ref. Data, 15, 1011–1071, 1986.
  - Kercher, J. P., Riedel, T. P., and Thornton, J. A.: Chlorine activation by N2O5: simultaneous, in situ detection of ClNO2 and N2O5 by chemical ionization mass spectrometry, Atmos. Meas. Tech., 2, 193–204, 2009.
- Koss, A. R., Warneke, C., Yuan, B., Coggon, M. M., Veres, P. R., and de Gouw, J. A.: Evaluation of NO+ reagent ion chemistry for online measurements of atmospheric volatile organic compounds, Atmos. Meas. Tech., 9, 2909–2925, https://doi.org/10.5194/amt-9-2909-2016, 2016.

- Koss, A. R., Sekimoto, K., Gilman, J. B., Selimovic, V., Coggon, M. M., Zarzana, K. J., Yuan, B., Lerner, B. M., Brown, S.
- S., Jimenez, J. L., Krechmer, J., Roberts, J. M., Warneke, C., Yokelson, R. J., and de Gouw, J.: Non-methane organic gas
- emissions from biomass burning: identification, quantification, and emission factors from PTR-ToF during the FIREX 2016

laboratory experiment, Atmos. Chem. Phys., 18, 3299–3319, https://doi.org/10.5194/acp-18-3299-2018, 2018.

- Krechmer, J., Lopez-Hilfiker, F., Koss, A., Hutterli, M., Stoermer, C., Deming, B., Kimmel, J., Warneke, C., Holzinger, R.,
- Jayne, J., Worsnop, D., Fuhrer, K., Gonin, M., and de Gouw, J.: Evaluation of a New Reagent-Ion Source and Focusing Ion-
- Molecule Reactor for Use in Proton-Transfer-Reaction Mass Spectrometry, Anal. Chem., 90, 12011-12018,
- 380 https://doi.org/10.1021/acs.analchem.8b02641, 2018.
- Lee, B. H., Lopez-Hilfiker, F. D., Mohr, C., Kurtén, T., Worsnop, D. R., and Thornton, J. A.: An Iodide-Adduct High-
  - Resolution Time-of-Flight Chemical-Ionization Mass Spectrometer: Application to Atmospheric Inorganic and Organic
  - Compounds, Environ. Sci. Technol., 48, 6309-6317, https://doi.org/10.1021/es500362a, 2014.
  - Lee, B. H., Lopez-Hilfiker, F. D., Veres, P. R., McDuffie, E. E., Fibiger, D. L., Sparks, T. L., Ebben, C. J., Green, J. R.,
- 385 Schroder, J. C., Campuzano-Jost, P., Iyer, S., D'Ambro, E. L., Schobesberger, S., Brown, S. S., Wooldridge, P. J., Cohen, R.
- C., Fiddler, M. N., Bililign, S., Jimenez, J. L., Kurtén, T., Weinheimer, A. J., Jaegle, L., and Thornton, J. A.: Flight deployment
  - of a high-resolution time-of-flight chemical ionization mass spectrometer: Observations of reactive halogen and nitrogen oxide
  - species, J. Geophys. Res.-Atmos., 123, 7670-7686, https://doi.org/10.1029/2017JD028082, 2018.
  - Lopez-Hilfiker, F. D., Mohr, C., Ehn, M., Rubach, F., Kleist, E., Wildt, J., Mentel, Th. F., Lutz, A., Hallquist, M., Worsnop,
- 390 D., and Thornton, J. A.: A novel method for online analysis of gas and particle composition: description and evaluation of a
  - Filter Inlet for Gases and AEROsols (FIGAERO), Atmos. Meas. Tech., 7, 983–1001, https://doi.org/10.5194/amt-7-983-2014,
  - 2014.
    - Lopez-Hilfiker, F. D., Iyer, S., Mohr, C., Lee, B. H., D'Ambro, E. L., Kurtén, T., and Thornton, J. A.: Constraining the
  - sensitivity of iodide adduct chemical ionization mass spectrometry to multifunctional organic molecules using the collision
- limit and thermodynamic stability of iodide ion adducts, Atmos. Meas. Tech., 9, 1505–1512, https://doi.org/10.5194/amt-9-1512.
- 1505-2016, 2016.
  - Lopez-Hilfiker, F. D., Pospisilova, V., Huang, W., Kalberer, M., Mohr, C., Stefenelli, G., Thornton, J. A., Baltensperger, U.,
  - Prevot, A. S. H., and Slowik, J. G.: An extractive electrospray ionization time-of-flight mass spectrometer (EESI-TOF) for
  - online measurement of atmospheric aerosol particles, Atmos. Meas. Tech., 12, 4867–4886, https://doi.org/10.5194/amt-12-
- 400 4867-2019, 2019.
  - Mielke, L. H., Furgeson, A., and Osthoff, H. D.: Observation of ClNO2 in a Mid-Continental Urban Environment, Environ.
  - Sci. Technol., 45, 8889–8896, https://doi.org/10.1021/es201955u, 2011.
  - Neuman, J. A., Huey, L. G., Dissly, R. W., Fehsenfeld, F. C., Flocke, F., Holecek, J. C., Holloway, J. S., Hubler, G., Jakoubek,
  - R., Nicks Jr., D. K., Parrish, D. D., Ryerson, T. B., Sueper, D. T., and Weinheimer, A. J.: Fast-response airborne in situ
- 405 measurements of HNO3 during the Texas 2000 Air Quality Study, J. Geophys. Res.-Atmos., 107, 4436,
- https://doi.org/10.1029/2001JD001437, 2002.

- Neuman, J. A., Trainer, M., Brown, S. S., Min, K.-E., Nowak, J. B., Parrish, D. D., Peischl, J., Pollack, I. B., Roberts, J. M., Ryerson, T. B., and Veres, P. R.: HONO emission and production determined from airborne measurements over the Southeast U.S., J. Geophys. Res. Atmos., 121, 9237–9250, https://doi.org/10.1002/2016JD025197, 2016.
- Osthoff, H. D., Roberts, J. M., Ravishankara, A. R., Williams, E. J., Lerner, B. M., Sommariva, R., Bates, T. S., Coffman, D., Quinn, P. K., Dibb, J. E., Stark, H., Burkholder, J. B., Talukdar, R. K., Meagher, J., Fehsenfeld, F. C., and Brown, S. S.: High levels of nitryl chloride in the polluted subtropical marine boundary layer, Nature Geosci, 1, 324–328, https://doi.org/10.1038/ngeo177, 2008.
  - Palm, B. B., Peng, Q., Fredrickson, C. D., Lee, B. H., Garofalo, L. A., Pothier, M. A., Kreidenweis, S. M., Farmer, D. K.,
- Pokhrel, R. P., Shen, Y., Murphy, S. M., Permar, W., Hu, L., Campos, T. L., Hall, S. R., Ullmann, K., Zhang, X., Flocke, F., Fischer, E. V., and Thornton, J. A.: Quantification of organic aerosol and brown carbon evolution in fresh wildfire plumes, Proc Natl Acad Sci USA, 202012218, https://doi.org/10.1073/pnas.2012218117, 2020.
  - Riedel, T. P., Bertram, T. H., Crisp, T. A., Williams, E. J., Lerner, B. M., Vlasenko, A., Li, S.-M., Gilman, J., de Gouw, J., Bon, D. M., Wagner, N. L., Brown, S. S., and Thornton, J. A.: Nitryl Chloride and Molecular Chlorine in the Coastal Marine
- 420 Boundary Layer, Environ. Sci. Technol., 46, 10463–10470, https://doi.org/10.1021/es204632r, 2012.
  - Sanchez, J., Tanner, D. J., Chen, D., Huey, L. G., and Ng, N. L.: A new technique for the direct detection of HO2 radicals using bromidechemical ionization mass spectrometry (Br-CIMS): initial characterization, Atmos. Meas. Tech., 9, 3851–3861, https://doi.org/10.5194/amt-9-3851-2016, 2016.
  - Slusher, D. L., Huey, L. G., Tanner, D. J., Flocke, F. M., and Roberts, J. M.: A thermal dissociation--chemical ionization mass
- spectrometry (TD-CIMS) technique for the simultaneous measurement of peroxyacyl nitrates and dinitrogen pentoxide, J. Geophys. Res., 109, D19315, https://doi.org/10.1029/2004JD004670, 2004.
  - Thornton, J. A., Kercher, J. P., Riedel, T. P., Wagner, N. L., Cozic, J., Holloway, J. S., Dubé, W. P., Wolfe, G. M., Quinn, P. K., Middlebrook, A. M., Alexander, B., and Brown, S. S.: A large atomic chlorine source inferred from mid-continental reactive nitrogen chemistry, Nature, 464, 271–274, https://doi.org/10.1038/nature08905, 2010.
- Veres, P., Roberts, J. M., Warneke, C., Welsh-Bon, D., Zahniser, M., Herndon, S., Fall, R., and de Gouw, J.: Development of negative-ion proton-transfer chemical-ionization mass spectrometry (NI-PT-CIMS) for the measurement of gas-phase organic acids in the atmosphere, International Journal of Mass Spectrometry, 274, 48–55, https://doi.org/10.1016/j.ijms.2008.04.032, 2008.
  - Veres, P. R., Neuman, J. A., Bertram, T. H., Assaf, E., Wolfe, G. M., Williamson, C. J., Weinzierl, B., Tilmes, S., Thompson,
- C. R., Thames, A. B., Schroder, J. C., Saiz-Lopez, A., Rollins, A. W., Roberts, J. M., Price, D., Peischl, J., Nault, B. A., Møller, K. H., Miller, D. O., Meinardi, S., Li, Q., Lamarque, J.-F., Kupc, A., Kjaergaard, H. G., Kinnison, D., Jimenez, J. L., Jernigan, C. M., Hornbrook, R. S., Hills, A., Dollner, M., Day, D. A., Cuevas, C. A., Campuzano-Jost, P., Burkholder, J. B., Bui, T. P., Brune, W. H., Brown, S. S., Brock, C. A., Bourgeois, I., Blake, D. R., Apel, E. C., and Ryerson, T. B.: Correction for Veres et al., Global airborne sampling reveals a previously unobserved dimethyl sulfide oxidation mechanism in the marine atmosphere,
- 440 Proc Natl Acad Sci USA, 118, e2113268118, https://doi.org/10.1073/pnas.2113268118, 2021.

Villalta, P. W. and Howard, C. J.: Direct Kinetics Study of the CH<sub>3</sub> C(O)O<sub>2</sub> + NO Reaction Using Chemical Ionization Mass Spectrometry, J. Phys. Chem., 100, 13624–13628, https://doi.org/10.1021/jp9614153, 1996.

Zaytsev, A., Breitenlechner, M., Koss, A. R., Lim, C. Y., Rowe, J. C., Kroll, J. H., and Keutsch, F. N.: Using collision-induced dissociation to constrain sensitivity of ammonia chemical ionization mass spectrometry (NH4+ CIMS) to oxygenated volatile organic compounds, Atmos. Meas. Tech., 12, 1861–1870, https://doi.org/10.5194/amt-12-1861-2019, 2019.

445